

**MUDr. Samer Asad**

*Gynecology clinic Prague -  
Obstetrics, treatment of infertility,  
in vitro fertilization (IVF)*



**MUDr. Matus Baran, MBA**

*Lecturer*



## Jett Plasma Lift **Medical** in gynecology

- Suspicious findings on the suppository
- Condylomata acuminata (genital warts)
- Herpes genitalis
- Vaginal rejuvenation
- Sutures after delivery (scars)

## Technical design – Device + applicators

- Contactless 1mm tip
- Golden tip with insulation
- Vaginal rejuvenation probe
- Special gynecology probe



# Pre-treatment protocol

Before start of the first treatment:

- Examination with PCR and Vulture cultivation
- If is present mykoplasma pneumoniae, Chlamydia infection, Ureaplasma urealyticum, etc., we cannot continue with any treatment



# Treatment protocol

■ Indication: **Suspicious findings on the suppository LSIL, ASCUS including HPV positive virus**

- 3x application of Jett Plasma for 8 min., intensity 5, within 6 weeks
- Local anaesthesia – only ultrasound gel
- Desinfection – Betadine
- Antioxidants - Folic acid, Gynerol, etc.

## HPV positive virus



# Treatment protocol


■ Indication: **Suspicious findings on the suppository LSIL, ASCUS including HPV positive virus**

- Treatment without any pain or restrictions
- Effect of Jett Plasma - mucosal stimulation, blood circulation, epithelium recovery
- healing a suppository which is resistant to viruses

# Treatment protocol

Indications **Condylomata, Herpes genitalis**

- 3x application of Jett Plasma for 8 min., intensity 5, within 6 weeks
- Local anaesthesia – lydocain 10%
- Desinfection – Betadine
- Antioxidants - Folic acid, Gynerol, etc.

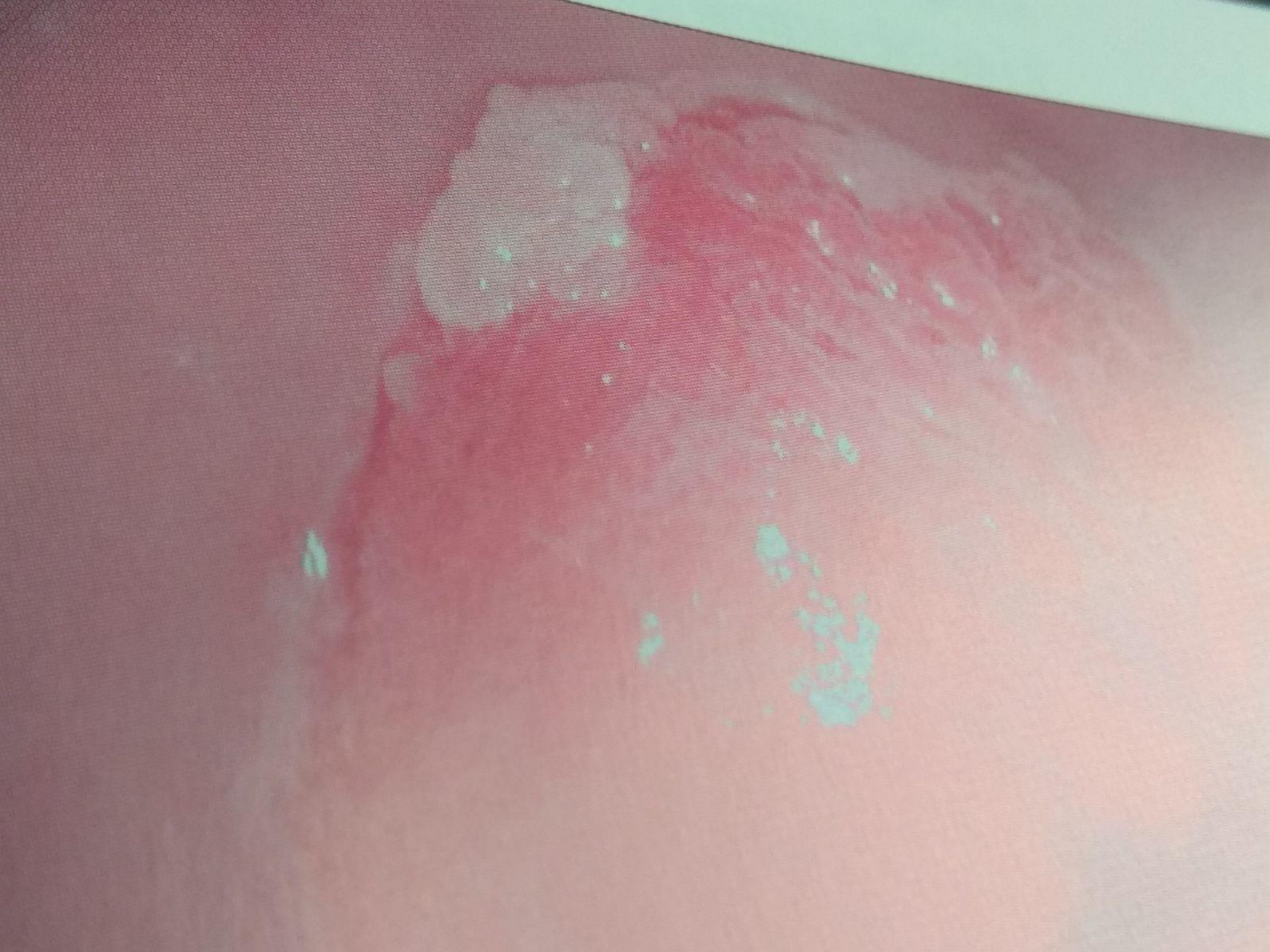

# Treatment protocol

■  
Indications **Condylomata, Herpes genitalis**

- Effect of Jett Plasma – ablation, removal of the bearings, necrotizing
- During and after the treatment is patient without any pain and restrictions

# Treatment protocol

Indication: **Vaginal rejuvenation**

- 3x application of Jett Plasma for 8 min., intensity 5, within 6 weeks
- Local anaesthesia – without
- Application of ultrasound gel

# Treatment protocol

■  
Indication: **Sutures after delivery (scars)**

- 3x application of Jett Plasma for 5-8 min., intensity 8, within 6 weeks
- Local anaesthesia – without
- Desinfection – Betadine (without alcohol)



# Experiences and feedbacks

- Condylomata – 20 patients
- Herpes genitalis – 11 patients
- Suspicious findings on the suppository –  
**56 patients**
- Vaginal rejuvenation – 10 patients
- Sutures after delivery (scars) – 2 patients

# Experiences and feedbacks

- All patients after the cytology control without finding
- All patients after the vaginal rejuvenation with better feelings
- Results without any restrictions and pain

# jett medical



---

■ [www.jett.eu](http://www.jett.eu) ■

Facebook: Jett Cosmetics & Jett Medical

Instagram: jett\_plasma\_lift\_medical