

PLASMA IN GYNECOLOGY

dermatology

GYNECOLOGY

ophthalmology



stomatology

ORL

proctology,
veterinary
medicine...



JETT PLASMA LIFT MEDICAL

- widespread use especially in dermatology
- unique method for blepharoplasty, removing of fibromas, haemangiomas, warts, scars etc.



TREATMENT

1. Nonisolated heads

the same indications as in dermatology

2. Isolated heads


- condylomata accuminata – caused by HPV 6,11
- HSV infections
- leukoplakia...


..deep in vagina



3. Gynecological heads

- vaginal rejuvenation
- incontinence

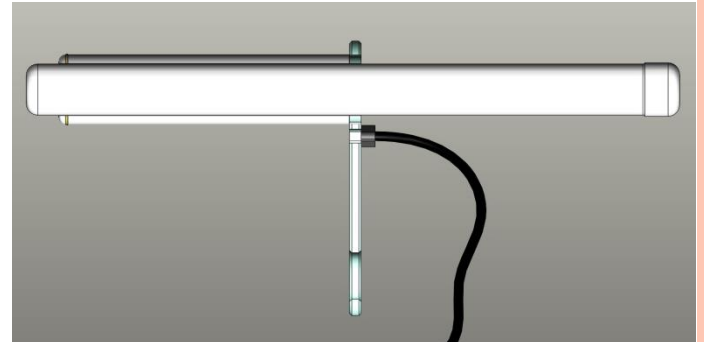

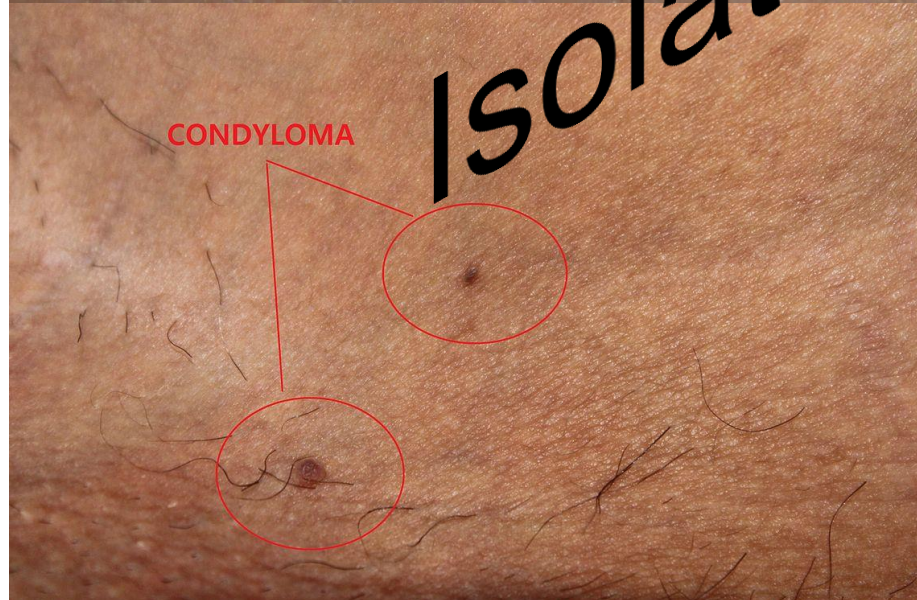



Isolated heads


GYNECOLOGICAL HEAD

- Our main procedure provided by plasma is VAGINAL REJUVENATION
- Vaginal Rejuvenation is widespread and popular method which solves problems with vaginal aging.
- It is currently performed as a surgery under general anaesthesia or as laser procedure.

Plasma should be safer than currently methods, comparison to laser plasma does not contain mirrors which can be moved in the device and cause adhesion or damage of tissue.



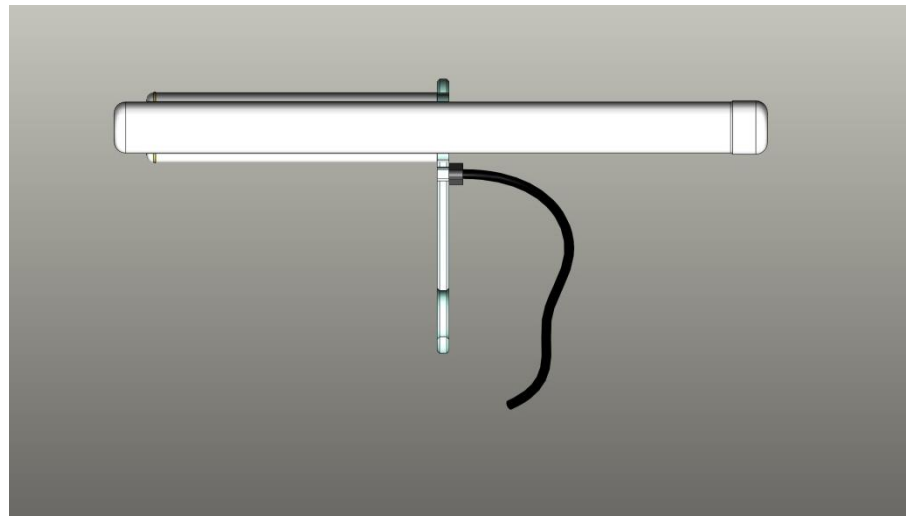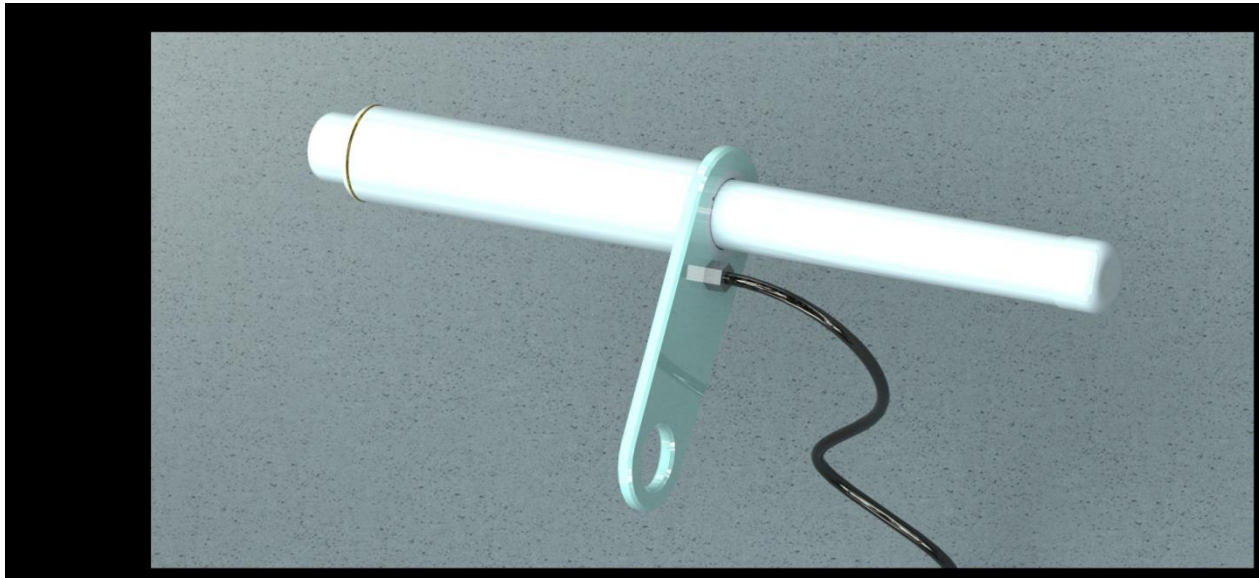

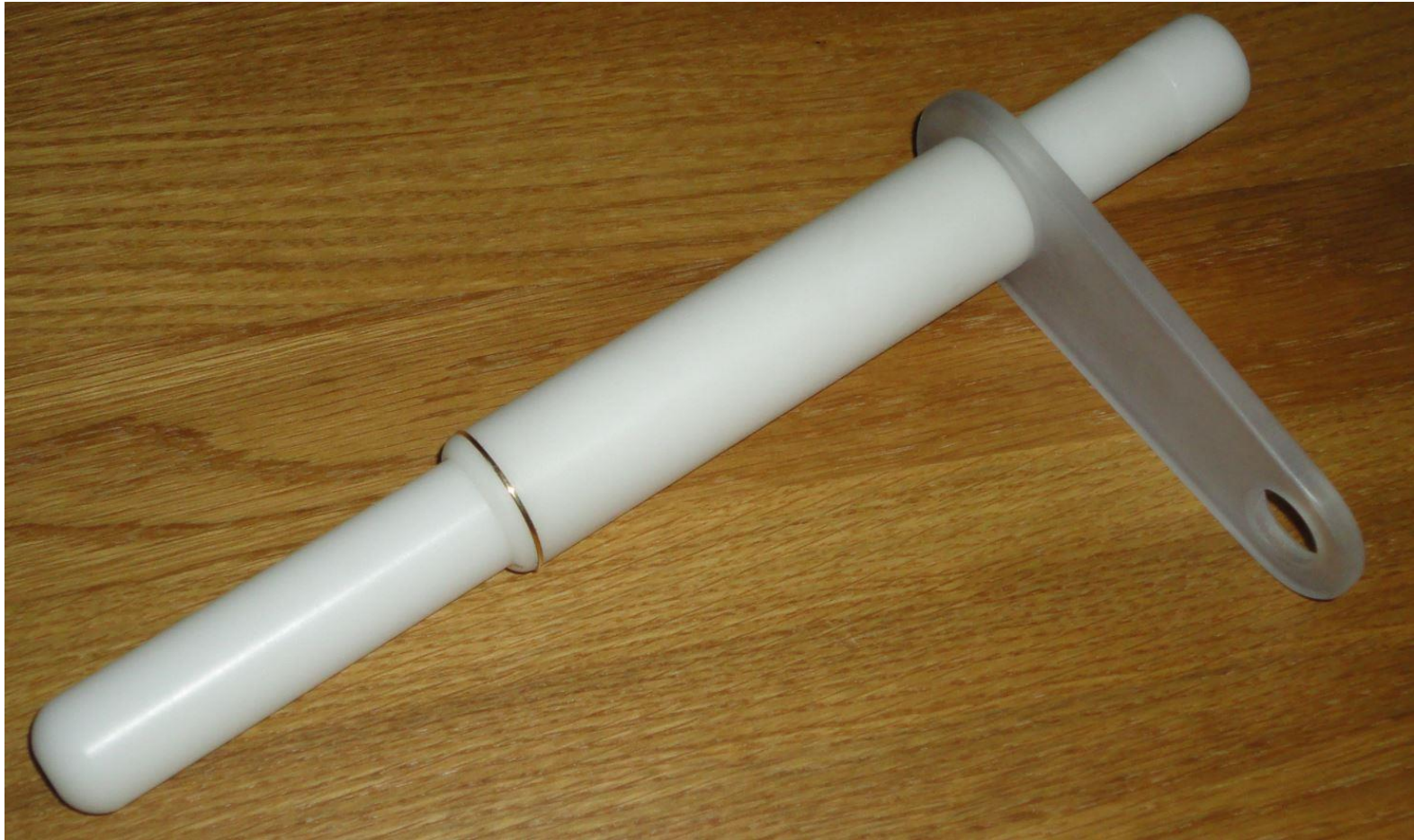

RECOMMENDED APPLICATION

- Now we are in stage of developing and trying it on patients, we cooperate with Czech and foreign doctors Because of disinfection, it is necessary to treat only one patient a day
- Treating includes:
 1. Insert the head into the vagina
 2. Set intensity 3 on device
 3. After this, turn the device on
 4. After this, slowly pull it out from the vagina – the recommended speed is 1 millimeter within 3seconds



- The head contains plasma source opened in one third of its circuit, the recommended intensity is 3
- Our goal is to have plasma source all opened circuit and work with higher intensity – 8
- If we will know exact information from doctors, we produce pneumatic system, which will make vaginal rejuvenation automatically



Thank you for your attention

